


Instructions Manual

Beauty instrument

WKB014

VACUUM & SPRAY

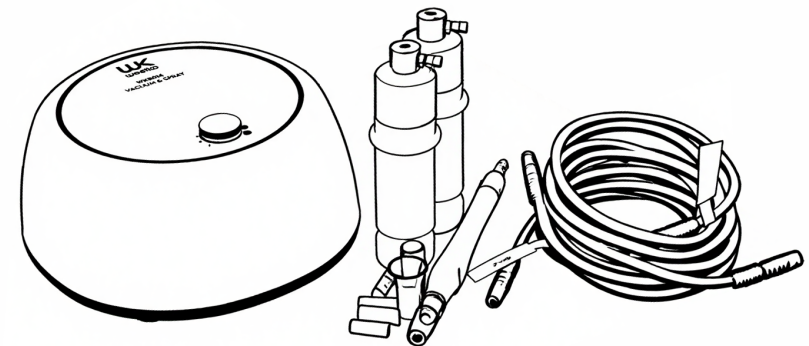

WK
Imported by
WEEKO BARCELONA SL

Pol. Ind. La Vall d'an
Camí de Garreta, 92
08600 Berga
Barcelona - Spain
T. +34 93 824 90 20
E. info@weelko.com

VERSION 1.0

Thank you sincerely for choosing our product. This manual is compiled to help users operate this instrument effectively. We hope our company's products and services will bring you a distinguished experience.

This user manual contains all general information about the instrument. It serves as the starting point for new users to begin using the instrument, as it provides a general introduction to the system, explains how it works, and details the instrument's installation and various operation methods.

Usage Precautions

1. Before turning on the machine, be sure to read the user manual carefully. Operators must be familiar with the user manual and operation methods, and the instrument must be managed by a designated person;
2. Do not use a power strip to avoid fire hazards;
3. When replacing the fuse, always turn off the power switch to prevent fire and electric shock;
4. Must use a power supply of 220V±22V, 50HZ±1HZ and an effectively grounded socket;
5. Using a damaged or modified power cord may cause fire and electric shock hazards;
6. If the instrument is loose, falling, or damaged, please contact the company in a timely manner;
7. Do not use the instrument in unstable environments such as tilting, vibrating, or impacting conditions;
8. No water or debris is allowed to enter the interior of the instrument. In case of foreign objects entering, please contact the company promptly;
9. When moving the instrument, the power plug must be pulled out first;
10. When moving the instrument, try to avoid vibration;
11. The instrument should be placed on a stable surface, with at least 5cm of space left around it to ensure air circulation and facilitate heat dissipation of the instrument;
12. Do not cover the instrument with anything to prevent blockage of the air vents;
13. Avoid using the instrument in the following environments: direct sunlight; high humidity environments; near water; vibrating or tilting places; strong magnetic field environments; places where electromagnetic waves or impulse voltages occur; places where chemicals are stored or corrosive gases exist;



Product Packaging List

Name	Image Quantity
Main Unit	1
Filter	1
Suction Pen (with 5 suction heads)	1 box
Spray Bottle	2
Power Cord	1

Product Specifications

- Input Voltage: AC 220V-240V
- Product Frequency: 50/60HZ
- Product Power: 20W
- Packaging Dimensions: 27x22x15CM
- Gross Weight: 1.6KG
- Net Weight: 1.35KG



1. Connect the device to a 220V power supply and press the power button. The device will emit a "beep" sound, indicating that it is powered on and enters standby mode;
2. For Suction Mode: Insert the suction tube into the filter nozzle, connect the suction pen to the suction tube, and insert the selected dermabrasion head. Press the small hole  on the handle to use it normally;
3. For Spray Mode: Insert the suction tube into the spray bottle connector and press the small hole of the air vent above the spray bottle to atomize the product;
4.  Suction Adjustment: Adjust the suction power from low to high using the knob.

14. This instrument shall not be placed or used near radios, televisions, copiers, fax machines, etc.;
15. Do not use it together with other instruments (such as microwaves) to avoid mutual interference between instruments and cause misoperation;
16. If an abnormal situation occurs to the instrument or the user during operation, the operator shall immediately move the operation handle away from the user, press the stop switch, and stop using the instrument;
17. The instrument shall be correctly used by skilled operators. Personnel other than operators are not allowed to use the instrument;
18. After the operation, place the control handle back on the control handle holder;
19. After using the operation handle, clean the probe thoroughly for next use.

Precautions for Troubleshooting, Storage, and Inspection

1. Do not use the instrument if it malfunctions. In case of smoke, burning smell, or inability to operate, continued use may cause fire or electric shock. Immediately press the stop button, turn off the power, unplug the power plug, and contact the company or distributor immediately;
2. Personnel other than professional technicians are strictly prohibited from repairing or modifying the instrument. The interior of the instrument must never be opened;
3. The instrument and its components must be inspected regularly;
4. Do not use cleaning tools that may damage the instrument's protective coating, such as brushes, brooms, etc.;
5. If the instrument is not used for more than 2 weeks, it must be carefully inspected before reuse.

Contraindicated Groups (Do Not Use If You Belong to the Following Groups)

- Individuals with implanted cardiac pacemakers, defibrillators, or metal implants (such as steel plates, bone nails) in the body;
- Pregnant women, lactating women, and women during menstruation;

- Patients with epilepsy, hypertension, severe heart disease, diabetes, tumors or cancer, or thyroid diseases;
- Individuals with open wounds, skin damage, severe inflammation, eczema, dermatitis on the skin, or those with facial telangiectasia (spider veins) or rosacea;
- Those who have received botulinum toxin injections, filler injections, laser treatments, chemical peels, or plastic surgery recently (usually within 6 months);
- Individuals with skin breakdown, inflammation, severe acne, eczema, rosacea, or those with facial metal implants.

Working Principle and Operation of the Vacuum & Spray Beauty Instrument

Working Principle

The vacuum & spray beauty instrument combines two core functions: negative pressure suction (vacuum) and product spraying (spray), providing dual skincare effects of deep cleansing and nutrient replenishment.

1. Deep Cleansing (Suction Mode)

It mainly uses vacuum negative pressure technology combined with special suction heads to extract dirt from pores.

- Preparation: Be sure to thoroughly clean the face first. Apply a warm towel to open the pores;
- Operation Points:
 - Select the Right Suction Head: Choose a suction head of appropriate size according to the cleaning area (such as T-zone, nasal alae);
 - Set Low Intensity: Start from the lowest suction level, especially for first-time use;
 - Correct Technique: Ensure the suction head is fully attached to the skin and sweep quickly. Do not stay in the same area for more than 2-3 seconds to avoid skin damage;
 - Avoid Areas: Strictly avoid the eye area, damaged, or red and swollen skin;
 - Cleaning and Maintenance: Immediately empty and clean the water tank and suction head after use, and keep them dry.

2. Essence Introduction/Soothe (Essence/Spray Mode)

The product adopts micro-hole atomization technology to convert liquid skincare products into fine mist particles, which evenly act on the surface of facial skin, replenishing moisture and nutrients for the skin and enhancing skin hydration. It is suitable for convenient skincare during daily home use and commuting breaks.

- Liquid Requirements: Use the special water tank of the instrument. Add purified water, toner, or special beauty essence. Do not add thick creams, essential oils, or irritating liquids containing alcohol;
- Water Filling Note: Do not exceed the water tank scale line when adding water to prevent liquid overflow and damage to internal parts of the instrument.

Core Safety Guidelines for Use

As a physically interfering instrument, improper operation of the vacuum & spray beauty instrument may damage the skin barrier. Please pay attention to the following:

- Frequency Control: Do not use it daily. It is recommended to use the deep cleansing mode 1-2 times a week; the specific frequency depends on individual skin type. For the introduction/care mode, refer to the product recommendations (e.g., the intensive care mode of LG BLP1 is recommended 3 times a week);
- Intensity and Duration: Prioritize lower intensity and shorter duration. Do not pursue high suction or prolonged stay on the same area;
- Post-Treatment Recovery Period: After facial plastic surgery, laser treatment, chemical peeling, and other projects, wait for the skin to fully recover (usually at least 2 months is recommended) before using the instrument;
- Special Reminder for Pregnant Women: Due to unstable hormone levels, pregnant or lactating women are not recommended to use the instrument.